


Ball Head Technology®

Pioneers in Angulation




PRODUCT CATALOGUE
CONVENTIONAL PROSTHESIS

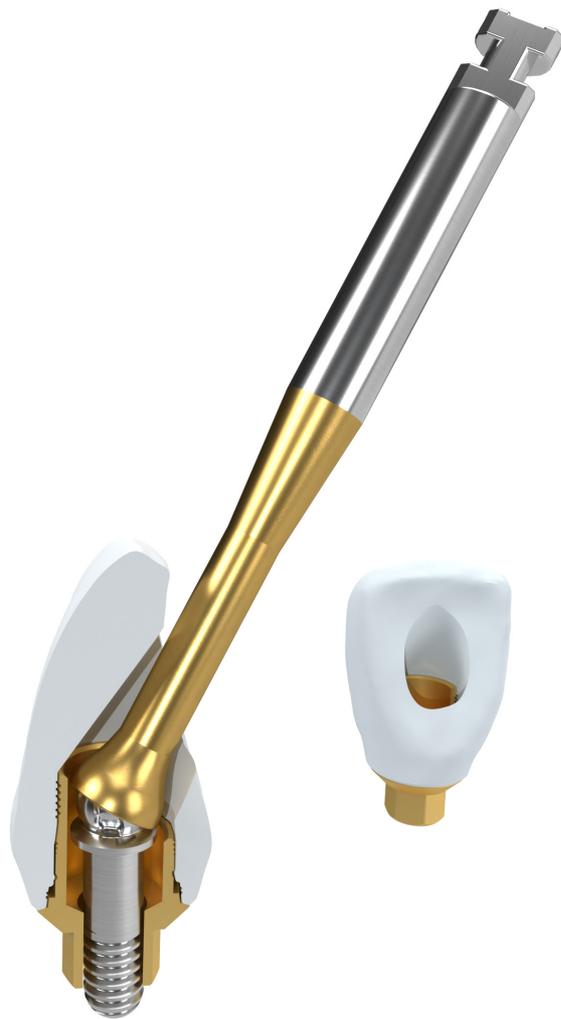

30° OF FREEDOM



The evolution towards customization in dental prosthetic design increasingly requires the use of angled abutments to provide versatility, simplicity and ease of use. Our Premium and Universal Solutions are pioneers in covering these needs.

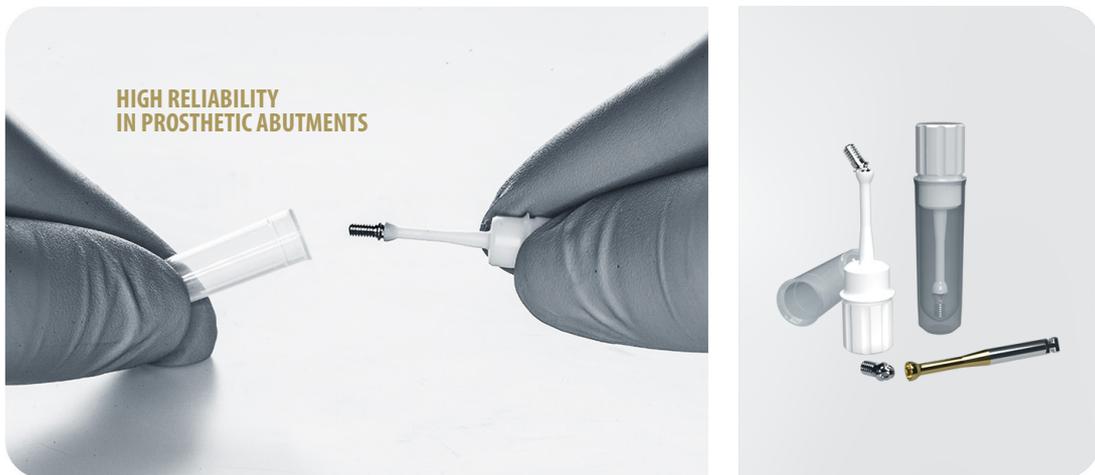

**HIGH RELIABILITY
IN PROSTHETIC ABUTMENTS**

Our main goal is achieving excellence to generate safety, reliability and comfort when using our products, both for the dental lab and clinical professionals, responsible for the patient's treatment.



**PIONEERS'
KNOW-HOW**

Since our first patent in 2008 we have developed the most efficient and innovative angled solutions based on the vision and experience of the most qualified restoration and oral surgery professionals.

INDEX

1.0 CONVENTIONAL PROSTHESIS	
1.1 DIRECT CASTING	5
1.2 OVERCASTING TO CrCo BASE	6
1.3 PRODUCTS	7/32
2.0 BHS30 SCREWDRIVER	33
3.0 PROTOCOLS	34/41
4.0 COMPATIBILITY CHART	42

Av. de Les Garrigues 17
25001 Lleida · Spain
(+34) 973 831 391
info@ballheadtechnology.com
www.ballheadtechnology.com

All names and trademarks used in this catalogue have been registered and are property of their respective companies.
Images published are merely intended for commercial purposes and have no contractual nature.
The information contained in this catalogue is aimed exclusively to professionals in the dental sector.
Ball Head Technology® reserves the right to make changes to this document without prior notice.
The guarantee conditions of Ball Head Technology® will apply as long as the instructions for each reference are respected.



1.1 DIRECT CASTING



SIMPLICITY EQUALS EFFICIENCY

- Engaging & non-engaging connections for all major implant systems.
- Continuous surface for an easy wax-up & casting.
- Optimized screw channel design and emergence profile.
- 10, 20 & 30° angle correction.
- Soft tissue care.

1.2 OVERCASTING TO CrCo BASE



PRECISION AND SAFETY FOR YOUR RESTORATIONS

- Continuous surface for an easy wax-up, casting & overcasting.
- Conical design for multiple restorations & engaging flat wall for single ones.
- Optimal emergence profile. Soft tissue care.
- Optimized access screw channel design.
- 10, 20 & 30° angle correction.
- Optimum height allowing 30° angulation for BHS screwdriver.
- Guide grip for BHS castable axial retention & free rotational movement.
- Retention pins & guide grip that enables axial retention.
- Engaging & non-engaging design for multiple and single restorations.

AB IMPLANTS COMPATIBLE

I SERIES INTERNAL HEX	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	NARROW Ø 3 - 3,75	30°	BHS07230	BHS07231	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	BHS08951A
		20°	BHS07232	BHS07233	BHS04057					
		10°	BHS07234	BHS07235	BHS04058					
	STANDARD Ø 3,5 - 6,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
10C CONUS 	RP Ø 4,2	30°	BHS07067	BHS07068	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	BHS08312A
		20°	BHS07069	BHS07070	BHS04105					
		10°	BHS07071	BHS07072	BHS04106					

ACE COMPATIBLE

INFINITY EXTERNAL HEX	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	NP Ø 3,3 / 3,75	30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP Ø 4,0	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064						
		10°	BHS07065	BHS07066						
INFINITY TRI-CAM 	NP Ø 3,5	30°	BHS07091	BHS07092				BHS02027	BHS02028	BHS08111A
		20°	BHS07093	BHS07094						
		10°	BHS07095	BHS07096						
	RP Ø 4,3	30°	BHS07016	BHS07017				BHS02009	BHS02010	BHS08112A
		20°	BHS07031	BHS07032						
		10°	BHS07046	BHS07047						
	WP Ø 5,0	30°	BHS07115	BHS07116				BHS02009	BHS02010	BHS08113A
		20°	BHS07117	BHS07118						
		10°	BHS07119	BHS07120						
INFINITY INTERNAL HEX 	NP Ø 3,7 / 4,1	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
	RP Ø 4,7 / 5,1	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039					
		10°	BHS07054	BHS07055	BHS04053					

ADIN COMPATIBLE

SWELL	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	Ø 3,5 - 6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
TOUAREG S/OS 	Ø 3,5 - 6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					

TOUAREG CLOSEFIT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT IMPLANTE
	Ø 3,5	30°		BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°		BHS07081	BHS07082	BHS04069					
		10°		BHS07083	BHS07084	BHS04070					
	Ø 4,3 / 5	30°		BHS07067	BHS07068	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	BHS08312A
		20°		BHS07069	BHS07070	BHS04105					
		10°		BHS07071	BHS07072	BHS04106					

ALPHABIO COMPATIBLE

EXTERNAL ALPHABIO		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°		BHS07012	BHS07013						
		10°		BHS07014	BHS07015						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WF Ø 5,1	30°		BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°		BHS07220	BHS07221						
		10°		BHS07222	BHS07223						
INTERNAL ALPHABIO		Ø 3,3 - 6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						

ALPHA DENT COMPATIBLE

ACTIVE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,75 / 4,2 / 5,0	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					
ACTIVE PLUS		Ø 4,2 / 5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS07073
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						
CLASSIC		Ø 3,3 / 3,75 / 4,2 / 5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						

ANCLADEN COMPATIBLE

ANCLALOCK		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,75 - 5	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					

ANTHOGRYR COMPATIBLE

ANTHOFIT HE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA	
		RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A	
			20°	BHS07004	BHS07005	BHS04005						
			10°	BHS07007	BHS07008	BHS04009						
		WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A	
			20°	BHS07063	BHS07064							
			10°	BHS07065	BHS07066							
		RP Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A	
			20°	BHS07035	BHS07036							
			10°	BHS07050	BHS07051							
		SYN_RP Ø 4,8	30°	BHS07022	BHS07023				BHS02013	BHS02014	BHS08713A	
			20°	BHS07037	BHS07038							
			10°	BHS07052	BHS07053							

ARDS COMPATIBLE

SMART		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		Ø 3,75 - 4,5	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					
		Ø 3,3 - 6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					
		Ø 3,3 - 6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					

ASTRA COMPATIBLE

AQUA		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		RP Ø 3.5 / 4.0	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
			20°	BHS07099	BHS07100	BHS04087					
			10°	BHS07101	BHS07102	BHS04088					
		WP Ø 4.5 / 5.0	30°	BHS07130	BHS07131				BHS02009	BHS02010	BHS08123A
			20°	BHS07132	BHS07133						
			10°	BHS07134	BHS07135						
		Ø 3,5	30°		BHS07103				BHS02003	BHS02004	BHS08124A
			20°		BHS07104						
			10°		BHS07105						
		Ø 3,5	30°		BHS07106				BHS02003	BHS02004	BHS08125A
			20°		BHS07107						
			10°		BHS07108						

EV	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	Ø 3,6	30°	BHS07269	BHS07270	BHS04056	BHS03154	BHS03155	BHS02053	BHS02054	BHS08152A
		20°	BHS07271	BHS07272	BHS04057					
		10°	BHS07273	BHS07274	BHS04058					
	Ø 4,2	30°	BHS07263	BHS07264	BHS04056	BHS03158	BHS03159	BHS02055	BHS02056	BHS08153A
		20°	BHS07265	BHS07266	BHS04057					
		10°	BHS07267	BHS07268	BHS04058					

AVINENT COMPATIBLE

HE	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	NP Ø 3,5	30°	BHS07010	BHS07011			BHS02031	BHS02032	BHS08611A	
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°	BHS07218	BHS07219			BHS02039	BHS02040	BHS08213A	
		20°	BHS07220	BHS07221						
		10°	BHS07222	BHS07223						
IH	Ø 3,5 / 4,1	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
MULTIUNIT	Ø 3,5	30°		BHS07003	BHS04002		BHS03008	BHS02003	BHS02004	BHS08411A
		20°		BHS07006	BHS04006					
		10°		BHS07009	BHS04010					
OCEAN CONICAL (C.C.)	Ø 3,5 Implant Ø 3,5 / Ø 4,0	30°	BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°	BHS07081	BHS07082	BHS04069					
		10°	BHS07083	BHS07084	BHS04070					
	Ø 4,1 Implant Ø 4,5 / Ø 5,0	30°	BHS07067	BHS07068	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	BHS08312A
		20°	BHS07069	BHS07070	BHS04105					
		10°	BHS07071	BHS07072	BHS04106					

BEGO COMPATIBLE

EH	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	NP Ø 2,7 / 3,1	30°	BHS07148	BHS07149			BHS02013	BHS02014	BHS08421A	
		20°	BHS07150	BHS07151						
		10°	BHS07152	BHS07153						
IH	Ø 3,75 / Ø 4,1	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					

BIOGENESIS COMPATIBLE

3ICON	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA	
	NP Ø 3,5	30°	BHS07194	BHS07195				BHS02001	BHS02002	BHS08611A	
		20°	BHS07196	BHS07197							
		10°	BHS07198	BHS07199							
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A	
		20°	BHS07004	BHS07005	BHS04005						
		10°	BHS07007	BHS07008	BHS04009						
	WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A	
		20°	BHS07063	BHS07064							
		10°	BHS07065	BHS07066							
	RP Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A	
		20°	BHS07035	BHS07036							
		10°	BHS07050	BHS07051							
	RP Ø 4,8	30°	BHS07022	BHS07023				BHS02013	BHS02014	BHS08713A	
		20°	BHS07037	BHS07038							
		10°	BHS07052	BHS07053							
	RP Ø 3,5 / 4,0	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A	
		20°	BHS07099	BHS07100	BHS04087						
		10°	BHS07101	BHS07102	BHS04088						
	WP Ø 4.5 / 5.0	30°	BHS07130	BHS07131				BHS02009	BHS02010	BHS08123A	
		20°	BHS07132	BHS07133							
		10°	BHS07134	BHS07135							
		RP Ø 3,5 / 4,0	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
			20°	BHS07099	BHS07100	BHS04087					
			10°	BHS07101	BHS07102	BHS04088					
WP Ø 4.5/ 5.0 / 5,5		30°	BHS07130	BHS07131				BHS02009	BHS02010	BHS08123A	
		20°	BHS07132	BHS07133							
		10°	BHS07134	BHS07135							
		RP / WP	30°		BHS07103				BHS02003	BHS02004	BHS08124A
			20°		BHS07104						
			10°		BHS07105						
	RP / WP	30°		BHS07106				BHS02003	BHS02004	BHS08125A	
		20°		BHS07107							
		10°		BHS07108							

BIOHORIZONS COMPATIBLE

EXTERNAL		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,7	30°		BHS07236	BHS07237				BHS02031	BHS02032	
		20°		BHS07238	BHS07239						
		10°		BHS07240	BHS07241						
<hr/>											
INTERNAL											
	NP Ø 3,5	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					
	RP Ø 4,5	30°		BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039					
		10°		BHS07054	BHS07055	BHS04053					
	WP Ø 5,7	30°		BHS07166	BHS07167				BHS02021	BHS02022	BHS08323A
		20°		BHS07168	BHS07169						
		10°		BHS07170	BHS07171						

BIOMET 3i COMPATIBLE

EH OSSEOTITE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07194	BHS07195				BHS02001	BHS02002	BHS08611A
		20°		BHS07196	BHS07197						
		10°		BHS07198	BHS07199						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WP Ø 5,0	30°		BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°		BHS07063	BHS07064						
		10°		BHS07065	BHS07066						
<hr/>											
CERTAIN											
	NP Ø 3,4	30°		BHS07029	BHS07030				BHS02007	BHS02008	BHS08511A
		20°		BHS07044	BHS07045						
		10°		BHS07059	BHS07060						
	RP Ø 4,1	30°		BHS07018	BHS07019				BHS02007	BHS02008	BHS08512A
		20°		BHS07033	BHS07034						
		10°		BHS07048	BHS07049						
	WP Ø 5,0	30°		BHS07121	BHS07122				BHS02007	BHS02008	BHS08513A
		20°		BHS07123	BHS07124						
		10°		BHS07125	BHS07126						

BIONER COMPATIBLE

IKELT / BIKELT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WP Ø 5,0	30°		BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°		BHS07063	BHS07064						
		10°		BHS07065	BHS07066						
	WP Ø 5,1	30°		BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°		BHS07220	BHS07221						
		10°		BHS07222	BHS07223						

HIKELT	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,8	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
	RP Ø 4,7	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039					
		10°	BHS07054	BHS07055	BHS04053					
TOP DM 	Ø 3,5/4,0/5,0	30°	BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°	BHS07081	BHS07082	BHS04069					
		10°	BHS07083	BHS07084	BHS04070					
MULTIUNIT 	RP Ø 4,8	30°		BHS07003	BHS04002		BHS03008	BHS02003	BHS02004	BHS08411A
		20°		BHS07006	BHS04006					
		10°		BHS07009	BHS04010					

BTI COMPATIBLE

HE	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Tiny Ø 3,5	30°	BHS07148	BHS07149				BHS02013	BHS02014	BHS08421A
		20°	BHS07150	BHS07151						
		10°	BHS07152	BHS07153						
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,5	30°	BHS07154	BHS07155				BHS02013	BHS02014	BHS08423A
		20°	BHS07156	BHS07157						
		10°	BHS07158	BHS07159						
CUATRILOBULAR 	RP Ø 4,1	30°	BHS07160	BHS07161				BHS02021	BHS02022	BHS08231A
		20°	BHS07162	BHS07163						
		10°	BHS07164	BHS07165						

BTK COMPATIBLE

8NECK INT	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	SR Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°	BHS07035	BHS07036						
		10°	BHS07050	BHS07051						
EVO DL EXT 	NP Ø 3,4/ Ø 3,7	30°	BHS07194	BHS07195				BHS02035	BHS02036	BHS08611A
		20°	BHS07196	BHS07197						
		10°	BHS07198	BHS07199						
	RP Ø 4,1/ Ø 4,6	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,0-5,6	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064						
		10°	BHS07065	BHS07066						

EVO DL INT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		Ø3,5 / Ø3,7 / Ø4,6 / Ø5,0 / Ø5,6	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					
KLASSIC EXT		NP Ø 3,4	30°	BHS07194	BHS07195				BHS02035	BHS02036	BHS08611A
	20°		BHS07196	BHS07197							
	10°		BHS07198	BHS07199							
		RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
			20°	BHS07004	BHS07005	BHS04005					
			10°	BHS07007	BHS07008	BHS04009					
		WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
			20°	BHS07063	BHS07064						
			10°	BHS07065	BHS07066						
KLASSIC INT		Ø3,5 / Ø4,5 / Ø5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						
KONIC EXT		NP Ø 3,4	30°	BHS07194	BHS07195				BHS02035	BHS02036	BHS08611A
	20°		BHS07196	BHS07197							
	10°		BHS07198	BHS07199							
		RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
			20°	BHS07004	BHS07005	BHS04005					
			10°	BHS07007	BHS07008	BHS04009					
		WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
			20°	BHS07063	BHS07064						
			10°	BHS07065	BHS07066						
KONIC INT		Ø 3,5 / 5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						
SAFE TL INT		SR Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
	20°		BHS07035	BHS07036							
	10°		BHS07050	BHS07051							

B&W COMPATIBLE

HEXAGONO EXT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		NP Ø 3,3	30°	BHS07136	BHS07137				BHS02035	BHS02036	BHS08741A
			20°	BHS07138	BHS07139						
			10°	BHS07140	BHS07141						
		RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
			20°	BHS07004	BHS07005	BHS04005					
			10°	BHS07007	BHS07008	BHS04009					
		WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
			20°	BHS07063	BHS07064						
			10°	BHS07065	BHS07066						
HEXAGONO INT		NP Ø 3,3	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
	20°		BHS07075	BHS07076	BHS04057						
	10°		BHS07077	BHS07078	BHS04058						
		RP Ø 4,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					
		WP Ø 5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					

CAMLOG COMPATIBLE

CAMLOG	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	Ø 3,3	30°	BHS07275	BHS07276	BHS04056	BHS03130	BHS03131	BHS02049	BHS02050	BHS08721A
		20°	BHS07277	BHS07278	BHS04057					
		10°	BHS07279	BHS07280	BHS04058					
	Ø 3,8	30°	BHS07281	BHS07282	BHS04056	BHS03134	BHS03135	BHS02049	BHS02050	BHS08722A
		20°	BHS07283	BHS07284	BHS04057					
		10°	BHS07285	BHS07286	BHS04058					
	Ø 4,3	30°	BHS07287	BHS07288	BHS04026	BHS03138	BHS03139	BHS02049	BHS02050	BHS08723A
		20°	BHS07289	BHS07290	BHS04039					
		10°	BHS07291	BHS07292	BHS04053					
	Ø 5,0	30°	BHS07293	BHS07294	BHS04065	BHS03142	BHS03143	BHS02051	BHS02052	BHS08724A
		20°	BHS07295	BHS07296	BHS04066					
		10°	BHS07297	BHS07298	BHS04067					

CONEXÃO SIST. PROTESE COMPATIBLE

MORSE INDEXADO	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	Ø3,5 / 3,75 / 4,0 / 5,0	30°	BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°	BHS07081	BHS07082	BHS04069					
		10°	BHS07083	BHS07084	BHS04070					
HEXAGONO EXT	Ø 3,3	30°	BHS07136	BHS07137				BHS02035	BHS02036	BHS08741A
		20°	BHS07138	BHS07139						
		10°	BHS07140	BHS07141						
	Ø 3,75/4,0	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064						
		10°	BHS07065	BHS07066						
HEXAGONO INT	Ø3,7 / 4,0/5,0	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039					
		10°	BHS07054	BHS07055	BHS04053					

CORTEX COMPATIBLE

DYNAMIX	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	Ø 3,3 / 3,8 / 4,2 / 5	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
CLASSIX	Ø 3,3 / 3,8 / 4,2 / 5	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
SATURN	Ø 3,8 / 4,2	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					

COWELLMEDI COMPATIBLE

EXTERNAL		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A	
		20°	BHS07004	BHS07005	BHS04005						
		10°	BHS07007	BHS07008	BHS04009						
	WP Ø 5,1	30°	BHS07218	BHS07219							
		20°	BHS07220	BHS07221							
		10°	BHS07222	BHS07223							

DENTALTECH COMPATIBLE

IMPLIC GII		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,5	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A	
		20°	BHS07039	BHS07040	BHS04039						
		10°	BHS07054	BHS07055	BHS04053						

IMPLASSIC HX2		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A	
		20°	BHS07004	BHS07005	BHS04005						
		10°	BHS07007	BHS07008	BHS04009						

DENTIS COMPATIBLE

S-CLEAN		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 3,4	30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A	
		20°	BHS07144	BHS07145	BHS04057						
		10°	BHS07146	BHS07147	BHS04058						

I-CLEAN		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A	
		20°	BHS07035	BHS07036							
		10°	BHS07050	BHS07051							

DENTIUM COMPATIBLE

SIMPLE LINE II		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RN Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A	
		20°	BHS07035	BHS07036							
		10°	BHS07050	BHS07051							

IMPLANTUM		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A	
		20°	BHS07144	BHS07145	BHS04057						
		10°	BHS07146	BHS07147	BHS04058						

SUPER LINE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		Ø 3,6 / 4 / 4,5 / 5 / 6	30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A
			20°	BHS07144	BHS07145	BHS04057					
			10°	BHS07146	BHS07147	BHS04058					

DENTSPLY SIRONA COMPATIBLE

XIVE S		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	D 3.4		30°	BHS07308	BHS07309	BHS04056	BHS03118	BHS03119	BHS02047	BHS02048	BHS08221
			20°	BHS07310	BHS07311	BHS04057					
			10°	BHS07312	BHS07313	BHS04058					
	D 3.8		30°	BHS07314	BHS07315	BHS04056	BHS03122	BHS03123	BHS02047	BHS02048	BHS08222A
			20°	BHS07316	BHS07317	BHS04057					
			10°	BHS07318	BHS07319	BHS04058					
	D 4.5		30°	BHS07320	BHS07321	BHS04056	BHS03126	BHS03127	BHS02047	BHS02048	BHS08223A
			20°	BHS07322	BHS07323	BHS04057					
			10°	BHS07324	BHS07325	BHS04058					

DIO IMPLANTS COMPATIBLE

EXTERNAL		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5		30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
			20°	BHS07012	BHS07013						
			10°	BHS07014	BHS07015						
	RP° Ø 4,1		30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
			20°	BHS07004	BHS07005	BHS04005					
			10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,1		30°	BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
			20°	BHS07220	BHS07221						
			10°	BHS07222	BHS07223						
UF II	Ø 3,8 / 4 / 4,5 / 5 / 5,5		30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A
			20°	BHS07144	BHS07145	BHS04057					
			10°	BHS07146	BHS07147	BHS04058					

EDISON MEDICAL COMPATIBLE

LAMINA (IH)II		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,3/3,75/4,2/5, 0/6,0		30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					
FORTIS (IH)	Ø 3,3/3,750 4,2/5,0/6,0		30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
			20°	BHS07075	BHS07076	BHS04057					
			10°	BHS07077	BHS07078	BHS04058					

ELITE MEDICA COMPATIBLE

EXTERNAL CONNECT	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	NP Ø 3,5	30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP ^o Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°	BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°	BHS07220	BHS07221						
		10°	BHS07222	BHS07223						


EUROTEKNIKA COMPATIBLE

UNEVA	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	Ø 3,6/4,1/4,8	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
NATURACTICS	Ø 3,5/4/4,5/5	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
NATURALL+	Ø 3,5/4/4,5/5	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
NATEA	NP Ø 3,6/4,1	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
	RP Ø 4,8/6	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
AESTHETICA	Ø 3,6/4,1/4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°	BHS07035	BHS07036						
		10°	BHS07050	BHS07051						

GALIMPLANT COMPATIBLE


IPXE - EXTERNA	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	RP Ø 3,5/4/5	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
ICI - INTERNA	RP Ø 3,2/3,5/4/5	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
IPX - INTERNA	RP Ø 3,5/4/4,5/5	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					

GMI COMPATIBLE

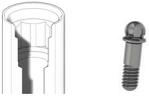FRONTIER (IH)	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	RP Standard Ø 3,3	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
	WP Wide Ø 4,3	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039					
		10°	BHS07054	BHS07055	BHS04053					
FRONTIER multi-aesthetic	RP Ø 4,8	30°	BHS07003	BHS07003	BHS04002	BHS03008	BHS02003	BHS02004	BHS08411A	
		20°	BHS07006	BHS07006	BHS04006					
		10°	BHS07009	BHS07009	BHS04010					
	WP Ø 6,0	30°	BHS07248	BHS07248	BHS04002	BHS02031	BHS02032	BHS02032	BHS02032	
		20°	BHS07249	BHS07249	BHS04006					
		10°	BHS07250	BHS07250	BHS04010					
PHOENIX (EH)	Narrow Ø 3,3	30°	BHS07236	BHS07237	BHS04001	BHS03003	BHS03004	BHS02031	BHS02032	BHS08612A
		20°	BHS07238	BHS07239						
		10°	BHS07240	BHS07241						
	Standard Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	Wide (Hex. 2,7) Ø 5,1	30°	BHS07061	BHS07062	BHS04005	BHS03003	BHS03004	BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064						
		10°	BHS07065	BHS07066						
	Wide (Hex. 3,4) Ø 5,1	30°	BHS07218	BHS07219	BHS04005	BHS03003	BHS03004	BHS02039	BHS02040	BHS08213A
		20°	BHS07220	BHS07221						
		10°	BHS07222	BHS07223						
	Extra-Wide Ø 6,0	30°	BHS07242	BHS07243	BHS04005	BHS03003	BHS03004	BHS02039	BHS02040	BHS08213A
		20°	BHS07244	BHS07245						
		10°	BHS07246	BHS07247						
FRONTIER multi-aesthetic	RP Ø 4,8	30°	BHS07003	BHS07003	BHS04002	BHS03008	BHS02003	BHS02004	BHS08411A	
		20°	BHS07006	BHS07006	BHS04006					
		10°	BHS07009	BHS07009	BHS04010					
	WP Ø 6,0	30°	BHS07248	BHS07248	BHS04002	BHS02031	BHS02032	BHS02032	BHS02032	
		20°	BHS07249	BHS07249	BHS04006					
		10°	BHS07250	BHS07250	BHS04010					

GT MEDICAL COMPATIBLE

BEST FIT HXE	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	NP Ø 3,5	30°	BHS07010	BHS07011	BHS04001	BHS03003	BHS03004	BHS02005	BHS02006	BHS08611A
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
WP Ø 5,1	30°	BHS07218	BHS07219	BHS04005	BHS03003	BHS03004	BHS02039	BHS02040	BHS08213A	
	20°	BHS07220	BHS07221							
	10°	BHS07222	BHS07223							
BEST FIT HXI	WP Ø 4,5 / 5,0	30°	BHS07130	BHS07131	BHS04005	BHS03003	BHS03004	BHS02009	BHS02010	BHS08123A
		20°	BHS07132	BHS07133						
		10°	BHS07134	BHS07135						

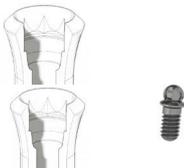BEST FIT 002		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					
	RN Ø 4,8	30°		BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°		BHS07035	BHS07036						
		10°		BHS07050	BHS07051						
	SYN_RN Ø 4,8	30°		BHS07022	BHS07023				BHS02013	BHS02014	BHS08713A
		20°		BHS07037	BHS07038						
		10°		BHS07052	BHS07053						

HI TECH COMPATIBLE

SELF THREAD		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	SLT Ø ,3 / 3,75	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					
	GTI Ø 4,2 / 5	30°		BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039					
		10°		BHS07054	BHS07055	BHS04053					
	Ø 3,5 / 4,3 / 5 / 6	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					

IMPLANT DIRECT COMPATIBLE

INTERACTIVE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°		BHS07081	BHS07082	BHS04069					
		10°		BHS07083	BHS07084	BHS04070					
	RP Ø 4,3	30°		BHS07067	BHS07068	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	BHS08312A
		20°		BHS07069	BHS07070	BHS04105					
		10°		BHS07071	BHS07072	BHS04106					
	NP Ø 3,5	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					
	RP Ø 4,5	30°		BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039					
		10°		BHS07054	BHS07055	BHS04053					
	WP Ø 5,7	30°		BHS07166	BHS07167				BHS02021	BHS02022	BHS08323A
		20°		BHS07168	BHS07169						
		10°		BHS07170	BHS07171						

REPLUS / REPLANT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07091	BHS07092				BHS02027	BHS02028	BHS08111A
		20°		BHS07093	BHS07094						
		10°		BHS07095	BHS07096						
	RP Ø 4,3	30°		BHS07016	BHS07017				BHS02009	BHS02010	BHS08112A
		20°		BHS07031	BHS07032						
		10°		BHS07046	BHS07047						
	WP Ø 5,0	30°		BHS07115	BHS07116				BHS02009	BHS02010	BHS08113A
		20°		BHS07117	BHS07118						
		10°		BHS07119	BHS07120						
SWISHPLUS / SWISHPLANT											
	RN Ø 4,8	30°		BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°		BHS07035	BHS07036						
		10°		BHS07050	BHS07051						
	SYN_RN Ø 4,8	30°		BHS07022	BHS07023				BHS02013	BHS02014	BHS08713A
		20°		BHS07037	BHS07038						
		10°		BHS07052	BHS07053						

DENTALCARE COMPATIBLE

JDEVOLUTION		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,7 / 4,3 / 5 / 6	30°		BHS07073	BHS07074	BHS04056			BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057	BHS03050	BHS03051			
		10°		BHS07077	BHS07078	BHS04058					


KEYSTONE DENTAL COMPATIBLE

RESTORE		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,3	30°		BHS07136	BHS07137				BHS02035	BHS02036	BHS08741A
		20°		BHS07138	BHS07139						
		10°		BHS07140	BHS07141						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001			BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005	BHS03003	BHS03004			
		10°		BHS07007	BHS07008	BHS04009					
OCT-XT											
	RN Ø 4,8	30°		BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°		BHS07035	BHS07036						
		10°		BHS07050	BHS07051						
	SYN_RN Ø 4,8	30°		BHS07022	BHS07023				BHS02013	BHS02014	BHS08713A
		20°		BHS07037	BHS07038						
		10°		BHS07052	BHS07053						
RENOVA											
	NP Ø 3,75	30°		BHS07073	BHS07074	BHS04056			BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057	BHS03050	BHS03051			
		10°		BHS07077	BHS07078	BHS04058					
	RP Ø 4,5 / 4,75	30°		BHS07024	BHS07025	BHS04026			BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039	BHS03033	BHS03034			
		10°		BHS07054	BHS07055	BHS04053					

KLOCKNER COMPATIBLE

SK2 (EH)	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,2	30° 20° 10°	BHS07109 BHS07111 BHS07113	BHS07110 BHS07112 BHS07114				BHS02031	BHS02032	BHS08521A
KL (EH)	NP Ø 3,5	30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001			BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005	BHS03003	BHS03004			
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°	BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°	BHS07220	BHS07221						
		10°	BHS07222	BHS07223						


LEADER COMPATIBLE

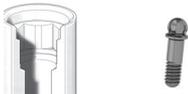TIXOS INTERNAL HEX	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5 / 4,0	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
TIXOS EXTERNAL HEX	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001			BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005	BHS03003	BHS03004			
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064						
		10°	BHS07065	BHS07066						

MEDENTIS COMPATIBLE

iCX	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,45 / 3,7 4,1 / 4,8	30° 20° 10°	BHS07257 BHS07259 BHS07261	BHS07258 BHS07260 BHS07262	BHS04056 BHS04057 BHS04058	BHS03180	BHS03181	BHS02015	BHS02016	BHS08751A

MEGAGEN COMPATIBLE

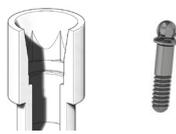
ANYONE INTERNAL		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA	
				⬡	◯		⬡	◯				
	NP		30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A	
			20°	BHS07144	BHS07145	BHS04057						
			10°	BHS07146	BHS07147	BHS04058						
	RP		30°	BHS07142	BHS07143	BHS04056						
			20°	BHS07144	BHS07145	BHS04057						
			10°	BHS07146	BHS07147	BHS04058						
	WP		30°	BHS07142	BHS07143	BHS04056						
			20°	BHS07144	BHS07145	BHS04057						
			10°	BHS07146	BHS07147	BHS04058						
ANYONE EXTERNAL												
	NP		30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A	
			20°	BHS07012	BHS07013							
			10°	BHS07014	BHS07015							
	RP		30°	BHS07001	BHS07002				BHS04001			
			20°	BHS07004	BHS07005				BHS04005			
			10°	BHS07007	BHS07008				BHS04009			
	WP		30°	BHS07061	BHS07062							
			20°	BHS07063	BHS07064							
			10°	BHS07065	BHS07066							
OCTA Abutment												
	SYN_RN Ø 4,8		30°	BHS07022	BHS07023	BHS04001	BHS03003	BHS03004		BHS02013	BHS02014	BHS08713A
			20°	BHS07037	BHS07038	BHS04005						
			10°	BHS07052	BHS07053	BHS04009						
ANYRIDGE												
	Small Ø 3,5		30°	BHS07251	BHS07252	BHS04056	BHS03166	BHS03167		BHS02027	BHS02028	BHS08991A
			20°	BHS07253	BHS07254	BHS04057						
			10°	BHS07255	BHS07256	BHS04058						
	Regular Ø 4,0/4,5		30°	BHS07251	BHS07252	BHS04056						
			20°	BHS07253	BHS07254	BHS04057						
			10°	BHS07255	BHS07256	BHS04058						
	Wide Ø 5,0/5,5		30°	BHS07251	BHS07252	BHS04056						
			20°	BHS07253	BHS07254	BHS04057						
			10°	BHS07255	BHS07256	BHS04058						

MICRODENT COMPATIBLE

EKTOS (IH)		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA	
				⬡	◯		⬡	◯				
	NP Ø 3,5		30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A	
			20°	BHS07075	BHS07076	BHS04057						
			10°	BHS07077	BHS07078	BHS04058						
	RP Ø 4,5		30°	BHS07024	BHS07025	BHS04026						
			20°	BHS07039	BHS07040	BHS04039						
			10°	BHS07054	BHS07055	BHS04053						
	WP Ø 5,7		30°	BHS07166	BHS07167							
			20°	BHS07168	BHS07169							
			10°	BHS07170	BHS07171							
MINI (EH)												
	NP Ø 3,5		30°	BHS07194	BHS07195							
			20°	BHS07196	BHS07197							
			10°	BHS07198	BHS07199							

SUV (EH)	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	SUV Ø4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A		
		20°	BHS07004	BHS07005	BHS04005							
		10°	BHS07007	BHS07008	BHS04009							
	SUV Ø5,1	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A		
		20°	BHS07063	BHS07064								
		10°	BHS07065	BHS07066								
	MST Ø3,5	30°	BHS07194	BHS07195				BHS02031	BHS02032	BHS08611A		
		20°	BHS07196	BHS07197								
		10°	BHS07198	BHS07199								
	MS Ø3,5	30°	BHS07194	BHS07195				BHS02031	BHS02032	BHS08611A		
		20°	BHS07196	BHS07197								
		10°	BHS07198	BHS07199								
	M/MV Ø4,2	30°	BHS07332	BHS07333				BHS02013	BHS02014	EM048014		
		20°	BHS07334	BHS07335								
		10°	BHS07336	BHS07337								
	MK Ø5,6	30°	BHS07242	BHS07243				BHS02039	BHS02040			
		20°	BHS07244	BHS07245								
		10°	BHS07246	BHS07247								
	MST Ø3,5	30°	BHS07194	BHS07195				BHS02031	BHS02032	BHS08611A		
		20°	BHS07196	BHS07197								
		10°	BHS07198	BHS07199								
	MS Ø3,5	30°	BHS07194	BHS07195				BHS02031	BHS02032	BHS08611A		
		20°	BHS07196	BHS07197								
		10°	BHS07198	BHS07199								
	SU Ø4,1	30°	BHS07001	BHS07002	BHS04001			BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005							
		10°	BHS07007	BHS07008	BHS04009							
	SU Ø5,1	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A		
		20°	BHS07063	BHS07064								
		10°	BHS07065	BHS07066								
	MK Ø5,6	30°	BHS07242	BHS07243				BHS02039	BHS02040			
		20°	BHS07244	BHS07245								
		10°	BHS07246	BHS07247								

MIS COMPATIBLE


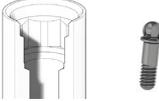
LANCE	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	STD Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A		
		20°	BHS07004	BHS07005	BHS04005							
		10°	BHS07007	BHS07008	BHS04009							
	WP Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A		
		20°	BHS07063	BHS07064								
		10°	BHS07065	BHS07066								
	NARROW Ø 3,30	30°	BHS07230	BHS07231	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	BHS08951A		
		20°	BHS07232	BHS07233	BHS04057							
		10°	BHS07234	BHS07235	BHS04058							
	STANDARD Ø 3,75 / 4,2	30°	BHS07073	BHS07074	BHS04056			BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057							
		10°	BHS07077	BHS07078	BHS04058							
	WIDE Ø 5,0 / 6,0	30°	BHS07024	BHS07025	BHS04026			BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039							
		10°	BHS07054	BHS07055	BHS04053							

MULTIUNIT MIS		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
		RP Ø 4,8	30° 20° 10°			BHS07127 BHS07128 BHS07129			BHS02037	BHS02038	BHS08411A

MOZO GRAU COMPATIBLE

INHEX		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
		RP Ø 3,75 / 4,25	30° 20° 10°			BHS07097 BHS07098 BHS07099 BHS07100 BHS07101 BHS07102	BHS04086 BHS04087 BHS04088	BHS03070 BHS03071	BHS02025	BHS02026	BHS08122A		
		WP Ø 5,0	30° 20° 10°			BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135		BHS02009	BHS02010	BHS08123A		
		CONO 20	30° 20° 10°				BHS07103 BHS07104 BHS07105		BHS02003	BHS02004	BHS08124A		
		CONO 45	30° 20° 10°				BHS07106 BHS07107 BHS07108		BHS02003	BHS02004	BHS08125A		
OSSEUS													
		NP Ø 3,4	30° 20° 10°				BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199		BHS02001	BHS02002	BHS08611A	
		RP Ø 4,1	30° 20° 10°			BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		WP Ø 5,0	30° 20° 10°				BHS07218 BHS07220 BHS07222	BHS07219 BHS07221 BHS07223		BHS02039	BHS02040	BHS08213A	

MPI COMPATIBLE

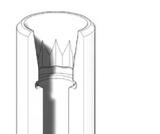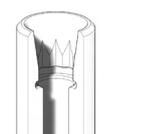EXTERNA EH		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
		NP Ø 3,5	30° 20° 10°						BHS02001	BHS02002	BHS08421A		
		RP Ø 4,1	30° 20° 10°			BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		WP Ø 5,0	30° 20° 10°				BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066		BHS02001	BHS02002	BHS08613A	
INTERNA CM													
		RP Ø 4,0	30° 20° 10°			BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		WP Ø 5,0	30° 20° 10°				BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135		BHS02009	BHS02010	BHS08123A	

NEOBIOTECH COMPATIBLE


EB EXTERNAL SYSTEM		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
				⬡	◯		⬡	◯			
	NP Ø 3,5	30°		BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°		BHS07012	BHS07013						
		10°		BHS07014	BHS07015						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°		BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°		BHS07220	BHS07221						
		10°		BHS07222	BHS07223						
IS INTERNAL SYSTEM											
	Ø 3,5 / 4 / 4,5 / 5	30°		BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A
		20°		BHS07144	BHS07145	BHS04057					
		10°		BHS07146	BHS07147	BHS04058					

NEODENT COMPATIBLE

TITAMAX HE SMART		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
				⬡	◯		⬡	◯			
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WP Ø 5,0	30°		BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A
		20°		BHS07063	BHS07064						
		10°		BHS07065	BHS07066						
MINI PILAR CM											
	RP Ø 4,8	30°			BHS07003	BHS04002		BHS03008	BHS02003	BHS02004	BHS08411A
		20°			BHS07006	BHS04006					
		10°			BHS07009	BHS04010					

NOBEL COMPATIBLE

BRÄNEMARK		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
				⬡	◯		⬡	◯			
	NP Ø 3,5	30°		BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°		BHS07012	BHS07013						
		10°		BHS07014	BHS07015						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°		BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°		BHS07220	BHS07221						
		10°		BHS07222	BHS07223						

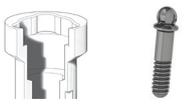REPLACE	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	NP Ø 3,5	30°	BHS07091	BHS07092						
		20°	BHS07093	BHS07094			BHS02027	BHS02028	BHS08111A	
		10°	BHS07095	BHS07096						
	RP Ø 4,3	30°	BHS07016	BHS07017						
		20°	BHS07031	BHS07032			BHS02009	BHS02010	BHS08112A	
		10°	BHS07046	BHS07047						
	WP Ø 5,0	30°	BHS07115	BHS07116						
		20°	BHS07117	BHS07118			BHS02009	BHS02010	BHS08113A	
		10°	BHS07119	BHS07120						
ACTIVE										
	NP Ø 3,5	30°	BHS07079	BHS07080	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	BHS08311A
		20°	BHS07081	BHS07082	BHS04069					
		10°	BHS07083	BHS07084	BHS04070					
	RP Ø 4,3 / Ø 5,0	30°	BHS07067	BHS07068	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	BHS08312A
		20°	BHS07069	BHS07070	BHS04105					
		10°	BHS07071	BHS07072	BHS04106					
MULTIUNIT	RP Ø 4,8	30°		BHS07003	BHS04002		BHS03008	BHS02003	BHS02004	BHS08411A
		20°		BHS07006	BHS04006					
		10°		BHS07009	BHS04010					

NT TRADING COMPATIBLE

2-CONnect	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA	
	M 1,4 Ø 5,0	30°		BHS07305	BHS04002	BHS03177	BHS02063	BHS02064	BHS08771A
		20°		BHS07306	BHS04006				
		10°		BHS07307	BHS04010				
	M 1,6 Ø 5,0	30°		BHS07305	BHS04002	BHS03177	BHS02065	BHS02066	BHS08771A
		20°		BHS07306	BHS04006				
		10°		BHS07307	BHS04010				
	M 1,8 Ø 5,0	30°		BHS07305	BHS04002	BHS03177	BHS02067	BHS02068	BHS08771A
		20°		BHS07306	BHS04006				
		10°		BHS07307	BHS04010				

OSSTEM IMPLANT COMPATIBLE

OSSTEM TS	PLATFORM Ø	∠	CASTABLE	CrCo CASTABLE	BASE CrCo	SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA		
	Ø 4,0 / 4,5 / 5,0 / 6,0	30°	BHS07142	BHS07143	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010	BHS08432A
		20°	BHS07144	BHS07145	BHS04057					
		10°	BHS07146	BHS07147	BHS04058					
OSSTEM US	MINI Ø 3,5	30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A
		20°	BHS07012	BHS07013						
		10°	BHS07014	BHS07015						
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
	WP Ø 5,1	30°	BHS07218	BHS07219				BHS02039	BHS02040	BHS08213A
		20°	BHS07220	BHS07221						
		10°	BHS07222	BHS07223						

OSTEOPLUS COMPATIBLE

ARGENTIUM		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
				⬡	◯		⬡	◯			
	RP Ø 4,2	30°	BHS07109	BHS07110					BHS02031	BHS02032	BHS08521A
		20°	BHS07111	BHS07112							
		10°	BHS07113	BHS07114							
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001		BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005						
		10°	BHS07007	BHS07008	BHS04009						
	WP Ø 5,0	30°	BHS07061	BHS07062					BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064							
		10°	BHS07065	BHS07066							
	RP Ø 3,5	30°	BHS07073	BHS07074	BHS04056		BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057						
		10°	BHS07077	BHS07078	BHS04058						
	WP Ø 4,5	30°	BHS07024	BHS07025	BHS04026		BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039						
		10°	BHS07054	BHS07055	BHS04053						

OXTEIN COMPATIBLE

HEX EXTERNA		PLATFOR M Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
				⬡	◯		⬡	◯			
	MINI Ø 3,3	30°	BHS07136	BHS07137					BHS02035	BHS02036	BHS08741A
		20°	BHS07138	BHS07139							
		10°	BHS07140	BHS07141							
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001		BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005						
		10°	BHS07007	BHS07008	BHS04009						
	WP Ø 5,0	30°	BHS07061	BHS07062					BHS02001	BHS02002	BHS08613A
		20°	BHS07063	BHS07064							
		10°	BHS07065	BHS07066							
	NP Ø 3,5	30°	BHS07073	BHS07074	BHS04056		BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057						
		10°	BHS07077	BHS07078	BHS04058						
	RP Ø 4,5	30°	BHS07024	BHS07025	BHS04026		BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039						
		10°	BHS07054	BHS07055	BHS04053						
	RN Ø 4,8	30°	BHS07020	BHS07021					BHS02011	BHS02012	BHS08712A
		20°	BHS07035	BHS07036							
		10°	BHS07050	BHS07051							
	RN Ø 4,8	30°	BHS07022	BHS07023					BHS02013	BHS02014	BHS08713A
		20°	BHS07037	BHS07038							
		10°	BHS07052	BHS07053							

PALTOP COMPATIBLE

INTERNAL HEX CONNECTION	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	NARROW Ø 3,25	30°	BHS07230	BHS07231	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	BHS08951A
		20°	BHS07232	BHS07233	BHS04057					
		10°	BHS07234	BHS07235	BHS04058					
	STANDARD Ø 3,75 / 4,2 / 5,0	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					
	WIDE Ø 6,0	30°	BHS07024	BHS07025	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A
		20°	BHS07039	BHS07040	BHS04039					
		10°	BHS07054	BHS07055	BHS04053					
PCA Implant System (Paltop Conical Active)										
	NARROW Ø 3,25	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
	STANDARD Ø 3,75/ 4,2 / 5,0	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A
		20°	BHS07099	BHS07100	BHS04087					
		10°	BHS07101	BHS07102	BHS04088					
MULTI UNIT										
	Ø 5,00	30°		BHS07326				BHS02031	BHS02032	BHS08271A
		20°		BHS07327						
		10°		BHS07328						
SINGLE UNIT										
	Ø 4,50	30°		BHS07329				BHS02031	BHS02032	BHS08251A
		20°		BHS07330						
		10°		BHS07331						

PHIBO COMPATIBLE

TSH	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	○		⬡	○			
	S2 Ø 3,3	30°	BHS07136	BHS07137				BHS02035	BHS02036	BHS08741A
		20°	BHS07138	BHS07139						
		10°	BHS07140	BHS07141						
	S3 / S4 Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
AUREA EVO										
	NP Ø 3,5	30°	BHS07178	BHS07179	BHS04056	BHS03106	BHS03107	BHS02041	BHS02042	
		20°	BHS07180	BHS07181	BHS04057					
		10°	BHS07182	BHS07183	BHS04058					
		0°			BHS04055					
	RP Ø 4,3	30°	BHS07186	BHS07187	BHS04056	BHS03110	BHS03111	BHS02021	BHS02022	
		20°	BHS07188	BHS07189	BHS04057					
		10°	BHS07190	BHS07191	BHS04058					
					BHS04055					
		0°								

RADHEX IMPLANTS COMPATIBLE

PHE / PHEA	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA						
			⬡	◯		⬡	◯									
	NP Ø 3,5	30°	BHS07010	BHS07011				BHS02005	BHS02006	BHS08611A						
		20°	BHS07012	BHS07013												
		10°	BHS07014	BHS07015												
	RP Ø 4,1	30°	BHS07001	BHS07002							BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005							BHS04005					
		10°	BHS07007	BHS07008							BHS04009					
PHI / PHIA	NP Ø 3,5	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A						
		20°	BHS07075	BHS07076	BHS04057											
		10°	BHS07077	BHS07078	BHS04058											
	RP Ø 4,5	30°	BHS07024	BHS07025	BHS04026						BHS03033	BHS03034	BHS02021	BHS02022	BHS08322A	
		20°	BHS07039	BHS07040	BHS04039											
		10°	BHS07054	BHS07055	BHS04053											
PCI	280	30°	BHS07097	BHS07098	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	BHS08122A						
		20°	BHS07099	BHS07100	BHS04087											
		10°	BHS07101	BHS07102	BHS04088											
	350	30°	BHS07130	BHS07131									BHS02009	BHS02010	BHS08123A	
		20°	BHS07132	BHS07133												
		10°	BHS07134	BHS07135												

SGS COMPATIBLE

P1 / P7	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	◯		⬡	◯			
	RP Ø 3,75	30°	BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°	BHS07075	BHS07076	BHS04057					
		10°	BHS07077	BHS07078	BHS04058					

SIGNO VINCES COMPATIBLE

DUO / INTEGRA	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
			⬡	◯		⬡	◯			
	RP Ø 4,1	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°	BHS07004	BHS07005	BHS04005					
		10°	BHS07007	BHS07008	BHS04009					
DUOCON / CM	WP_4.5/ 5.0	30°	BHS07130	BHS07131				BHS02009	BHS02010	BHS08123A
		20°	BHS07132	BHS07133						
		10°	BHS07134	BHS07135						

STRAUMANN COMPATIBLE

OCTA INTERN	PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA					
			⬡	◯		⬡	◯								
	RN Ø 4,8	30°	BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A					
		20°	BHS07035	BHS07036											
		10°	BHS07050	BHS07051											
SYN_RN Ø 4,8	30°	BHS07022	BHS07023										BHS02013	BHS02014	BHS08713A
	20°	BHS07037	BHS07038												
	10°	BHS07052	BHS07053												

BONE LEVEL		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NC Ø 3,3	30°		BHS07085	BHS07086				BHS02015	BHS02016	BHS08811A
		20°		BHS07087	BHS07088						
		10°		BHS07089	BHS07090						
	RC Ø 4,1 / 4,8	30°		BHS07027	BHS07028				BHS02015	BHS02016	BHS08812A
		20°		BHS07042	BHS07043						
		10°		BHS07057	BHS07058						
BL MULTI BASE		RC Ø 4,8	30°		BHS07026				BHS02017	BHS02018	
	20°			BHS07041							
	10°			BHS07056							

SWEDEN & MARTINA COMPATIBLE

PREMIUM KOHNO		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,3	30°		BHS07200	BHS07201				BHS02021	BHS02022	BHS08961A
		20°		BHS07202	BHS07203						
		10°		BHS07204	BHS07205						
	RP Ø 3,8	30°		BHS07206	BHS07207				BHS02021	BHS02022	BHS08962A
		20°		BHS07208	BHS07209						
		10°		BHS07210	BHS07211						
	WP Ø 4,25	30°		BHS07212	BHS07213				BHS02011	BHS02012	BHS08963A
		20°		BHS07214	BHS07215						
		10°		BHS07216	BHS07217						
SYRA (EH)		Ø 4,10 / 4,25	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
	20°		BHS07004	BHS07005	BHS04005						
	10°		BHS07007	BHS07008	BHS04009						
	Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A	
		20°	BHS07063	BHS07064							
		10°	BHS07065	BHS07066							
SYRA SHORT (EH)		Ø 4,10	30°	BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
	20°		BHS07004	BHS07005	BHS04005						
	10°		BHS07007	BHS07008	BHS04009						
	Ø 5,0	30°	BHS07061	BHS07062				BHS02001	BHS02002	BHS08613A	
		20°	BHS07063	BHS07064							
		10°	BHS07065	BHS07066							

TRI DENTAL IMPLANTS COMPATIBLE

TRI-VENT		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	Ø 3,75 / 4,1 / 4,7	30°		BHS07073	BHS07074	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057					
		10°		BHS07077	BHS07078	BHS04058					

TRINON COMPATIBLE

TRINON - Q2		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005					
		10°		BHS07007	BHS07008	BHS04009					

TRINON - QK		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	RN Ø 4,8	30°		BHS07020	BHS07021				BHS02011	BHS02012	BHS08712A
		20°		BHS07035	BHS07036						
		10°		BHS07050	BHS07051						
	NC Ø 3,3	30°		BHS07022	BHS07023			BHS02013	BHS02014	BHS08713A	
		20°		BHS07037	BHS07038						
		10°		BHS07052	BHS07053						

ZIACOM COMPATIBLE

ZINIC - IH		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,2	30°		BHS07230	BHS07231	BHS04056			BHS02041	BHS02042	BHS08951A
		20°		BHS07232	BHS07233	BHS04057	BHS03200	BHS03201			
		10°		BHS07234	BHS07235	BHS04058					
	RP Ø 3,5	30°		BHS07073	BHS07074	BHS04056			BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057	BHS03050	BHS03051			
		10°		BHS07077	BHS07078	BHS04058					
	WP Ø 4,5	30°		BHS07024	BHS07025	BHS04026			BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039	BHS03033	BHS03034			
		10°		BHS07054	BHS07055	BHS04053					
OEX - EH	NP Ø 3,3	30°		BHS07136	BHS07137				BHS02035	BHS02036	BHS08741A
		20°		BHS07138	BHS07139						
		10°		BHS07140	BHS07141						
	RP Ø 4,1	30°		BHS07001	BHS07002	BHS04001			BHS02001	BHS02002	BHS08612A
		20°		BHS07004	BHS07005	BHS04005	BHS03003	BHS03004			
		10°		BHS07007	BHS07008	BHS04009					
OIN - IH (M2)	RP Ø 4,0	30°		BHS07338	BHS07339	BHS04056			BHS02011	BHS02012	BHS08321A
		20°		BHS07340	BHS07341	BHS04057	BHS03050	BHS03051			
		10°		BHS07342	BHS07343	BHS04058					
	WP Ø 5,0	30°		BHS07338	BHS07339	BHS04026			BHS02011	BHS02012	BHS08322A
		20°		BHS07340	BHS07341	BHS04039	BHS03033	BHS03034			
		10°		BHS07342	BHS07343	BHS04053					
XDrive (Multi-unit)	RP Ø 4,8	30°		BHS07003	BHS07003	BHS04002			BHS02003	BHS02004	BHS08411A
		20°		BHS07006	BHS07006	BHS04006		BHS03008			
		10°		BHS07009	BHS07009	BHS04010					

ZIMMER COMPATIBLE

INTERNAL HEX		PLATFORM Ø	∠	CASTABLE		CrCo CASTABLE	BASE CrCo		SCREW + HANDLER	SCREW BHS	IMPLANT REPLICA
	NP Ø 3,5	30°		BHS07073	BHS07074	BHS04056			BHS02021	BHS02022	BHS08321A
		20°		BHS07075	BHS07076	BHS04057	BHS03050	BHS03051			
		10°		BHS07077	BHS07078	BHS04058					
	RP Ø 4,5	30°		BHS07024	BHS07025	BHS04026			BHS02021	BHS02022	BHS08322A
		20°		BHS07039	BHS07040	BHS04039	BHS03033	BHS03034			
		10°		BHS07054	BHS07055	BHS04053					
	WP Ø 5,7	30°		BHS07166	BHS07167				BHS02021	BHS02022	BHS08323A
		20°		BHS07168	BHS07169						
		10°		BHS07170	BHS07171						
SWISS-PLUS	RP Ø 4,8	30°		BHS07020	BHS07021				BHS02011	BHS02012	EM45144
		20°		BHS07035	BHS07036						
		10°		BHS07050	BHS07051						
	RP Ø 4,8	30°		BHS07022	BHS07023				BHS02013	BHS02014	EM46375
		20°		BHS07037	BHS07038						
		10°		BHS07052	BHS07053						

2.0 BHS30 SCREWDRIVER



EXCELLENT MECHANICAL PROPERTIES & SCREW HEAD PROTECTION

Universal contra angle compatible connection (ISO1797).

Hardened Steel offering excellent torsion resistance (90 N-cm).

Optimal coupling engagement ensuring screw head integrity.

TiN Coated: Surface treatment which increases the life of the driver.

Designed to correct 30° angle deviations with 100% efficient torque transmission.

SCREWDRIVER BHS



BHS01001 - L25mm
BHS01004 - L20mm
BHS01005 - L32mm

HANDLER



BHS01110 - Nob Ø 8mm / CA
BHS01111 - Str Ø 7mm / CA

3.0 PROTOCOLS


DIRECT CASTING

Products Required:


BHS30 Screw
BHS30 Screwdriver
Angled Castable Cylinder 10|20|30°

Materials:


Casting in CrCo, CrNi, Ti, Pd, Au,...

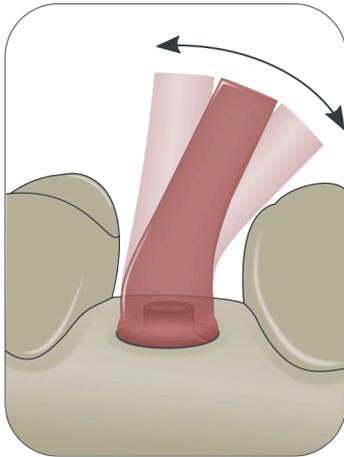


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration.

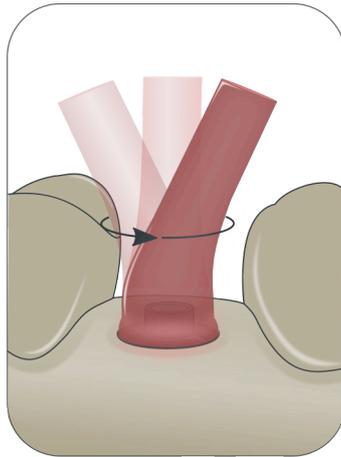


2.
LAB_ Insertion and orientation of the Angled Castable Cylinder on the analogue to obtain an angled access channel with optimum design and emergence profile. Fixation of the abutment with BHS30 Screw.

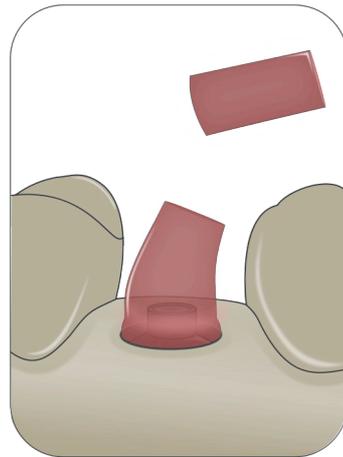


3.
LAB_ Cutting and adjustment of the Angled Castable Cylinder in accordance with the anatomy of the prosthesis and its antagonist.

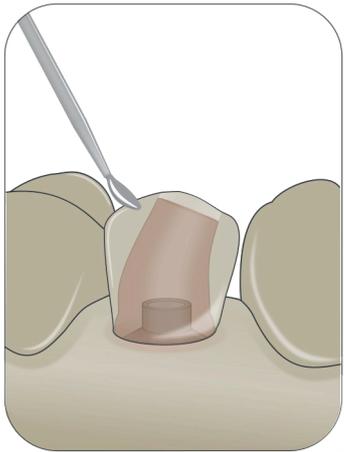

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Cylinder.

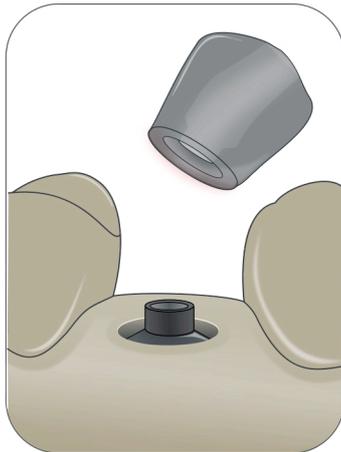


5.
LAB_ Metal casting of the structure and hand finishing of its implant connection to improve fitting on the implant analogue.

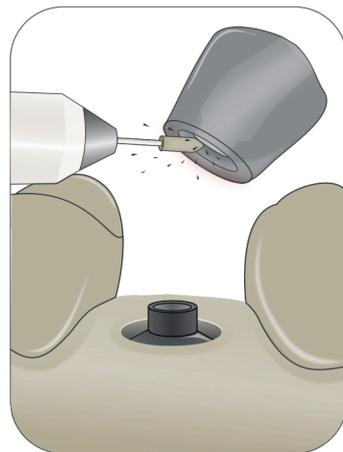


6.
LAB_ Metal casting of the structure and hand finishing of its implant connection to improve fitting on the implant analogue.

3.0 PROTOCOLS


DIRECT CASTING

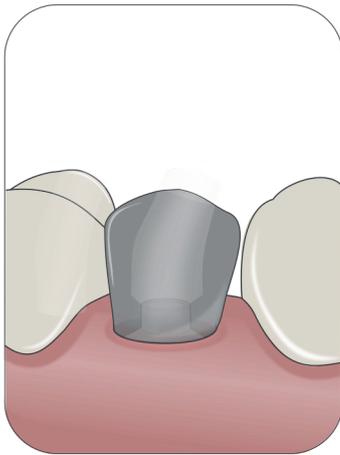

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth.

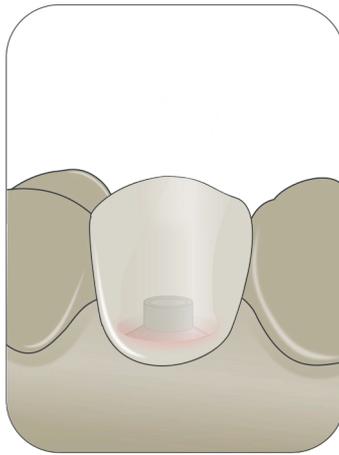

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.

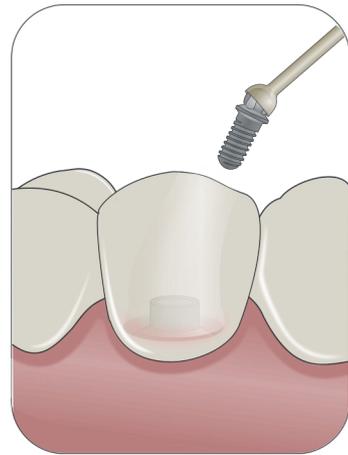

9.
CLINIC_ Placement of the prosthesis to the patient.


3.0 PROTOCOLS

OVERCASTING ON BHS30 CrCo BASE


Products Required:

BHS30 Screw
BHS30 Screwdriver
BHS30 CrCo Base
Angled Castable Chimney 10|20|30°


Materials:
Overcasting in CrCo

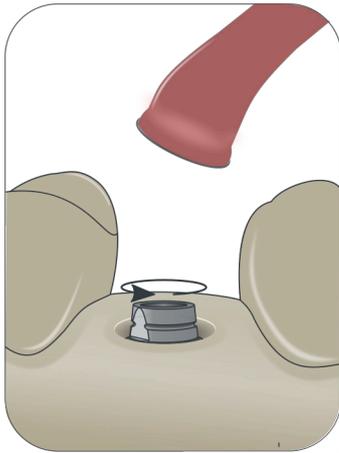


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 CrCo Base.

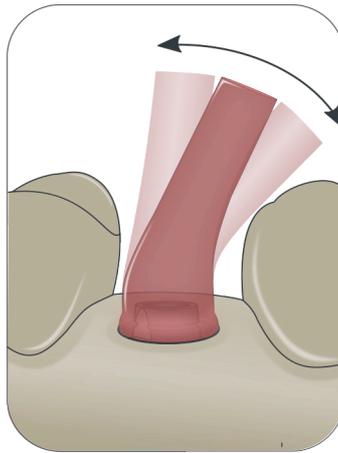


2.
LAB_ Insertion of the Base BHS30 CrCo + Chimney on the implant analogue to obtain an angled access channel with optimal design. Screwing of the BHS30 CrCo Base and final adjustment of the Chimney's orientation according to the desired position for screw access.

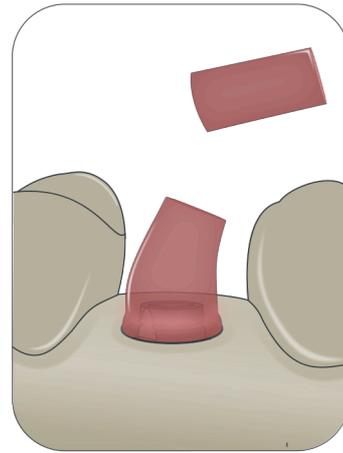


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney in accordance with the anatomy of the prosthesis and its antagonist.

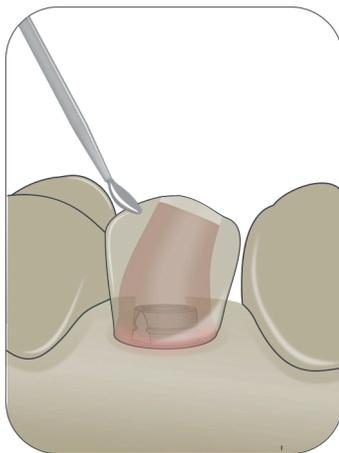

4.
LAB_ Wax-up and modelling of the structure on the BHS30 CrCo BASE + Angled Castable Chimney.

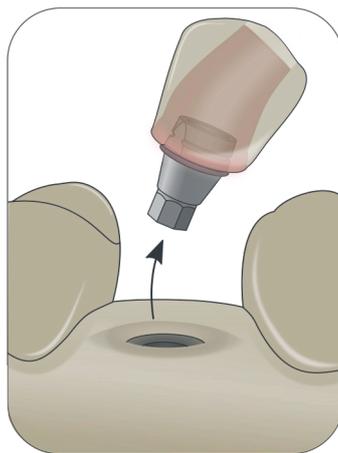


5.
LAB_ Metal overcasting of the structure.

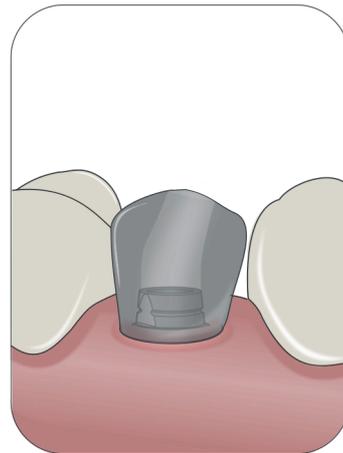


6.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth.

3.0 PROTOCOLS


OVERCASTING ON BHS30 CrCo BASE

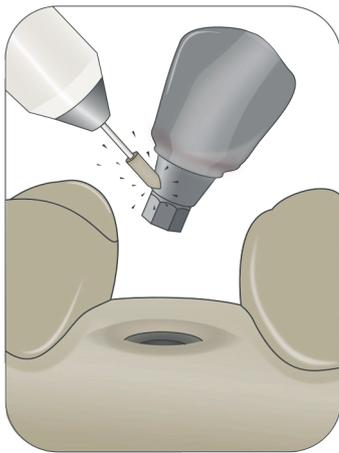

7.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.

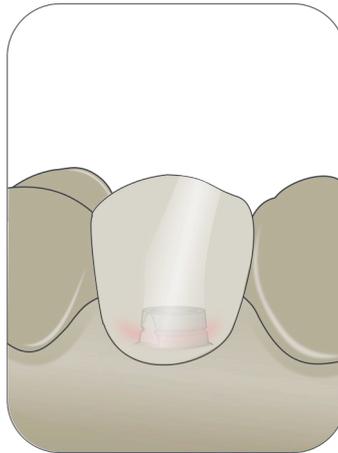

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.

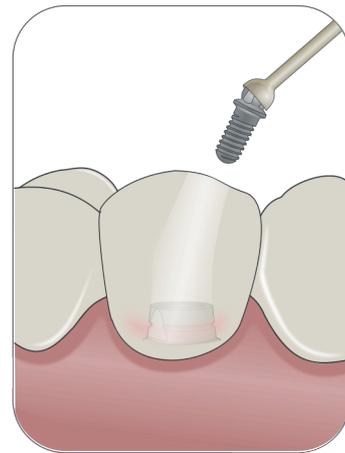

9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS


CASTING + WELDING ON BHS30 CrCo BASE

Products Required:


BHS30 Screw
BHS30 Screwdriver
BHS30 CrCo Base
Angled Castable Chimney 10|20|30°

Materials:


Casting in CrCo

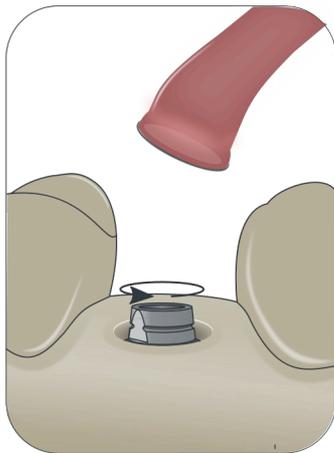


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 CrCo Base.

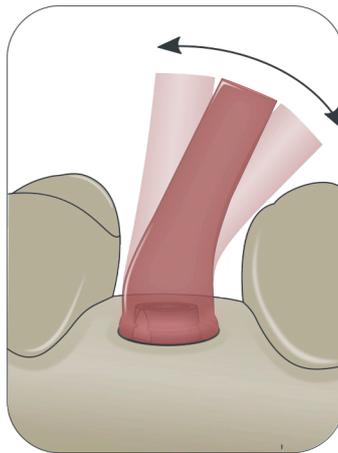


2.
LAB_ Insertion of the Base BHS30 CrCo + Chimney on the implant analogue to obtain an angled access channel with optimal design. Screwing of the BHS30 CrCo Base and final adjustment of the Chimney's orientation according to the desired position for screw access.

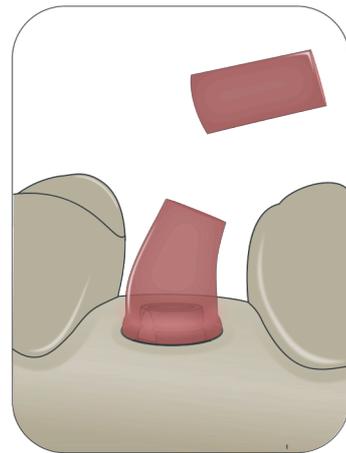


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney in accordance with the anatomy of the prosthesis and its antagonist.

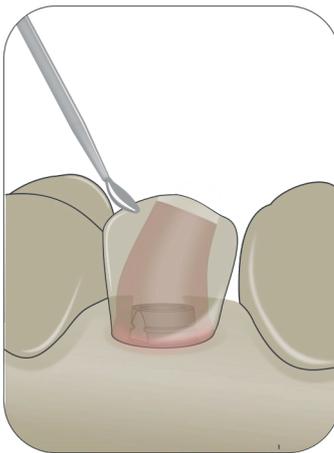

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Chimney.

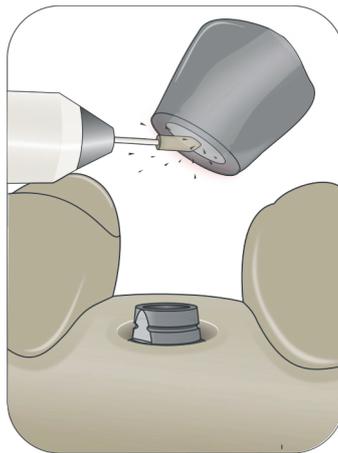


5.
LAB_ Removal of the Angled Castable Chimney with the modelling of the structure for metal casting, leaving the BHS30 CrCo Base screwed onto the model.

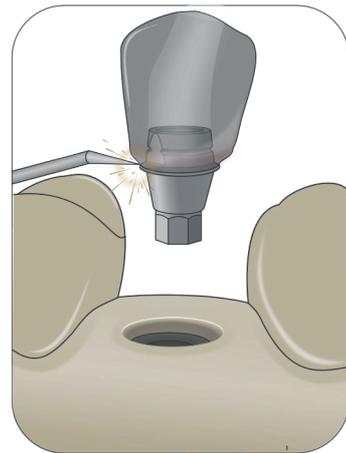


6.
A. LAB_ Adjustment of the metal cast structure on the BHS30 CrCo Base until obtaining passive fit on the model and perimetral laser welding of both parts and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS


CASTING + WELDING ON BHS30 CrCo BASE

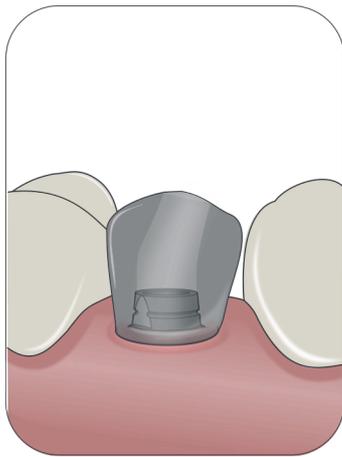

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.

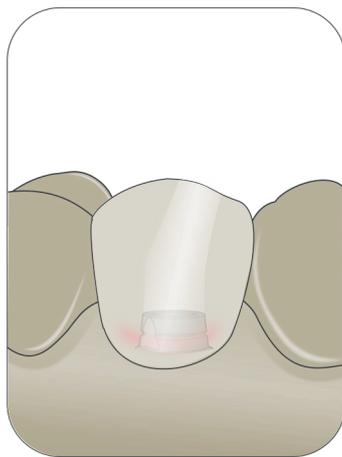

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and perimetral laser welding of both parts. Aesthetic restoration, polishing and finishing.

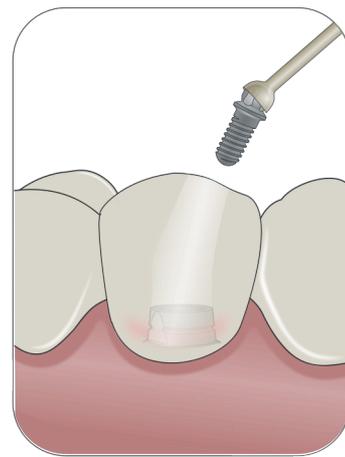

9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS


CASTING + CEMENTING ON BHS30 Ti BASE

Products Required:


BHS30 Screw
BHS30 Screwdriver
BHS30 Ti Base
Angled Castable Chimney 10|20|30°

Materials:


Casting in metal (Ti, CrCo, CrNi, noble...)

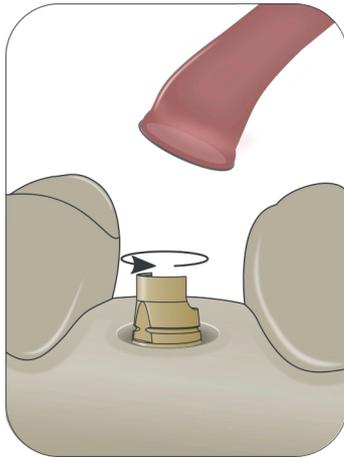


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 Ti Base.

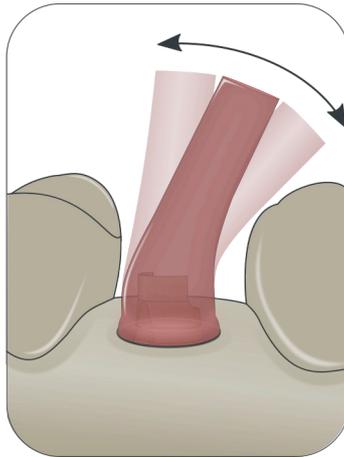


2.
LAB_ Insertion and orientation of the BHS30 Ti Base + Chimney on the implant analogue to obtain an angled access channel with optimal design and emergence profile. Screwing of the BHS30 Ti Base.

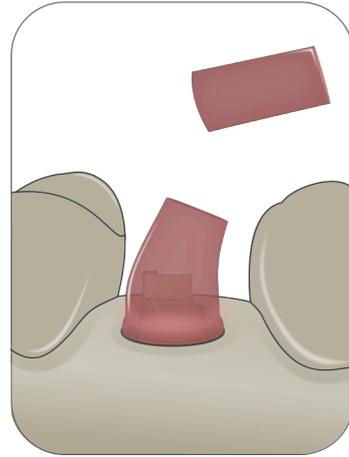


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney according to anatomy of the prosthesis and its antagonist.

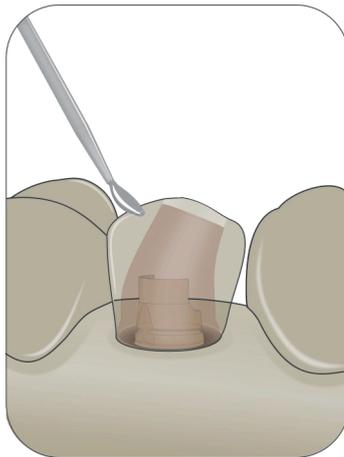

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Chimney.

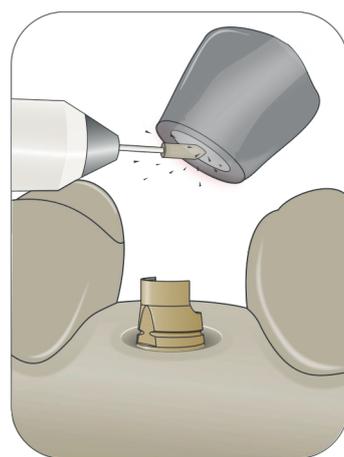


5.
LAB_ Removal of the Angled Castable Chimney with the modelling of the structure for metal casting, leaving the BHS30 Ti Base screwed onto the model.

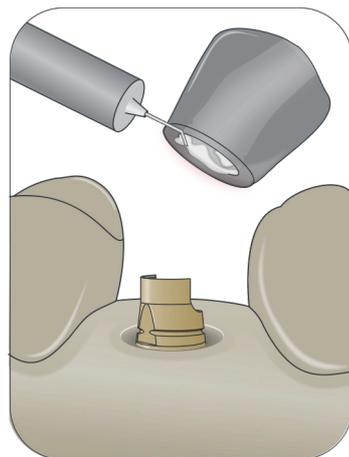


6.
A. LAB_ Adjustment of the metal cast structure on the BHS30 Ti Base until obtaining passive fit on the model and cementing of both parts, the excess of cement being expelled through the relining holes in the Chimney and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS


CASTING + CEMENTING ON BHS30 Ti BASE

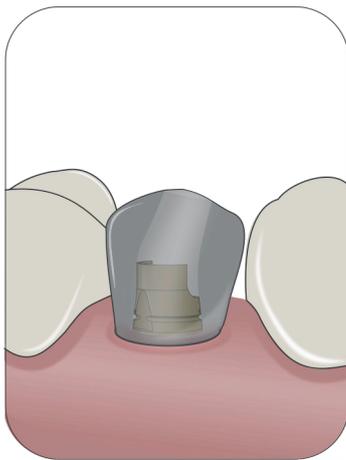

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.

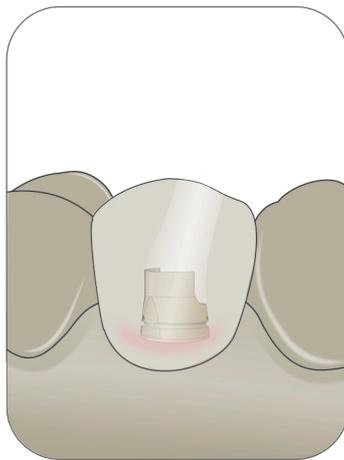

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and cementing of both parts, the excess of cement being expelled through the relining holes in the Chimney. Aesthetic restoration, polishing and finishing.

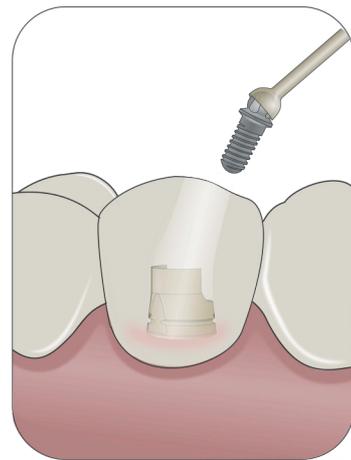

9.
CLINIC_ Placement of the prosthesis to the patient.

4.0 COMPATIBILITY CHART

TRADEMARK	PAGES	TRADEMARK	PAGES
AB IMPLANTS	7	IMPLANT DIRECT	20
ACE	7	JDENTAL CARE	21
ADIN	7	KEYSTONE DENTAL	21
ALPHABIO	8	KLOCKNER	22
ALPHADENT	8	LEADER	22
ANCLADEN	8	MEDENTIS	22
ANTHOGYR	9	MEGAGEN	23
ARDS	9	MICRODENT	23
ASTRA	9	MIS	24
AVINENT	10	MOZO GRAU	25
BEGO	10	MPI	25
BIOGENESIS	11	NEOBIOTECH	26
BIOHORIZONS	12	NEODENT	26
BIOMET 3i	12	NOBEL	26
BIONER	12	NT TRADING	27
BTI	13	OSSTEM IMPLANT	27
BTK	13	OSTEOPLUS	28
B&W	14	OXTEIN	28
CAMLOG	15	PALTOP	29
CON. SISTEMA PROTESE	15	PHIBO	29
CORTEX	15	RADHEX IMPLANTS	30
COWELLMEDI	16	SGS	30
DENTALTECH	16	SIGNO VINCES	30
DENTIS	16	STRAUMANN	30
DENTIUM	16	SWEDEN & MARTINA	31
DENTSPLY SIRONA	17	TRI DENTAL IMPLANTS	31
DIO IMPLANTS	17	TRINON	31
EDISON MEDICAL	17	ZIACOM	32
ELITE MEDICA	18	ZIMMER	32
EUROTEKNIKA	18		
GALIMPLANT	18		
GMI	19		
GT MEDICAL	19		
HI TECH	20		



www.ballheadtechnology.com

(+34) 973 831 391 · info@ballheadtechnology.com